



# Search My Commands Install Guide

---

L-3770

Last update: June 2014



## INTRODUCTION

This document covers the main installation steps for the Dragon 360 Medical | Search My Commands and other products if applicable (Clinical Calculators and PowerPack Content).

This document describes specifically the installation into the Dragon Network Edition of 2.1 or higher. It is intended to guide an NMC administrator through the process.

## PRE-REQUISITES

- Dragon Medical 360 | Network Edition 2.0 or higher.
- Administrator license to the Nuance Management Console and a member of the Administrator Group.

Features	Minimum	Recommended
RAM	2GB	4GB
Cache	512KB L2	2MB L2
Processor	Pentium IV or AMD Athlon 64-1GHz	
32-bit OS	XP Home/Professional SP2 or Higher; Vista Home/Pro SP1 or Higher; Window 7	
64-bit OS	XP Not supported; Vista Home/Pro SP1; Windows 7	
	Microsoft.NET Framework 4 (for v2 SMC)	

## DVD

The delivered DVD contains all files needed for installation of the products purchased. The DVD also contains documents pertaining to the products.

## iSupport

*iSupport Solution 16430: DM360 PowerPack™ / Search My Commands™ - Software Download and Documentation*


## Dragon 360 Medical | Search My Commands

### Installation Guide


(Within Dragon Network Edition 2.+)

#### Installation Steps:

1. Open the Nuance Management Console and log in.
2. Click on “Groups “ from the ribbon menu., and then “Add” under the groups section

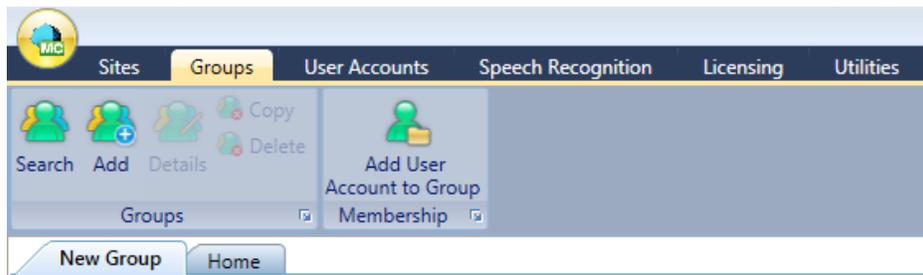

3. In the New Group Tab fill in the “Group Name”, “Site”, and “Role”

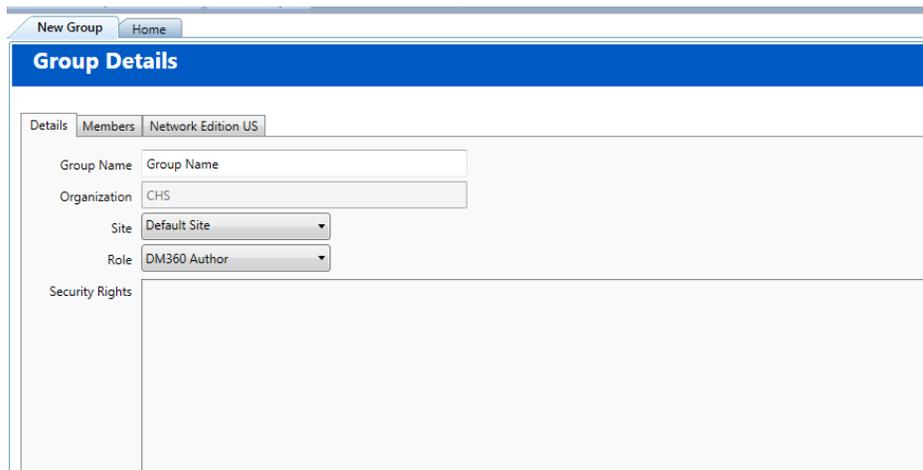


4. Click on “Members” tab and add at least one user to the group.
5. Click on “Save” to save the group.




- Click on “Speech Recognition” from ribbon menu.




- Click “Manage” within Command Sets section of the ribbon.
- Click on the correct group: Command Sets > Site > group within the main window (on the left).

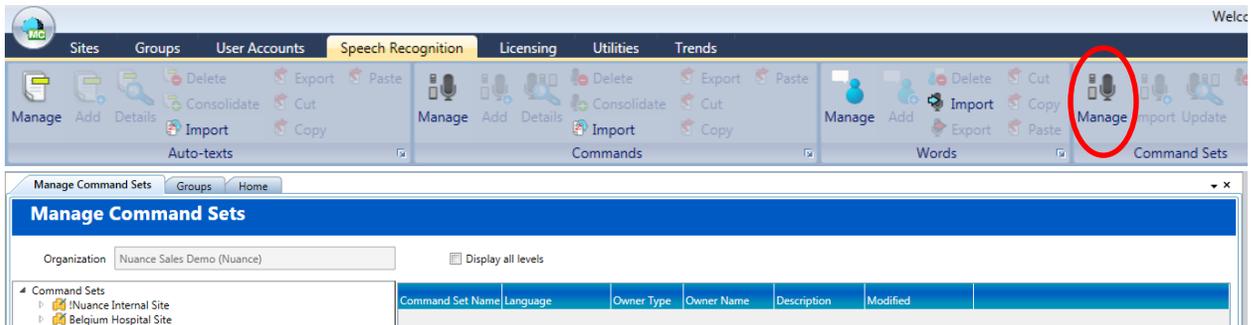


- From the Ribbon Menu click on “Import”

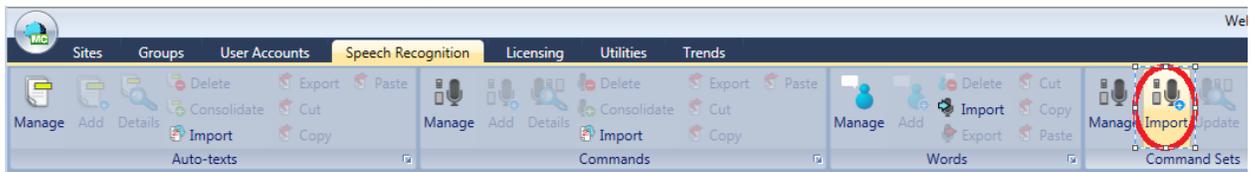


- In the “Import Command Set” Window, Click the Search Icon indicated with a magnifying glass.

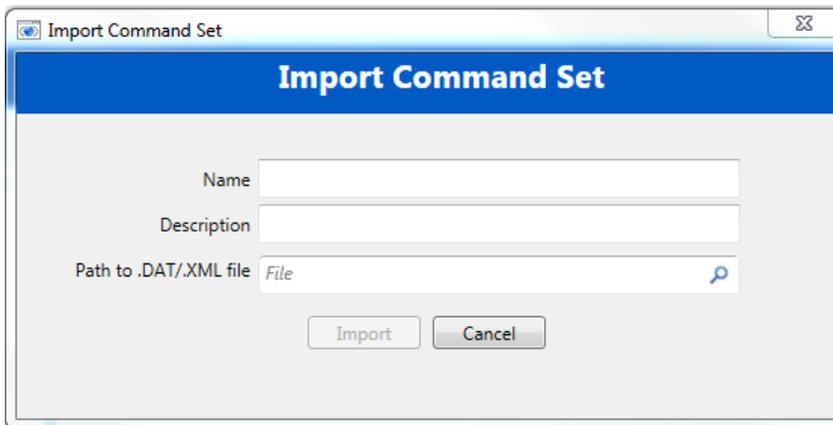

11. Select the Dat file to Import and Click “Open”

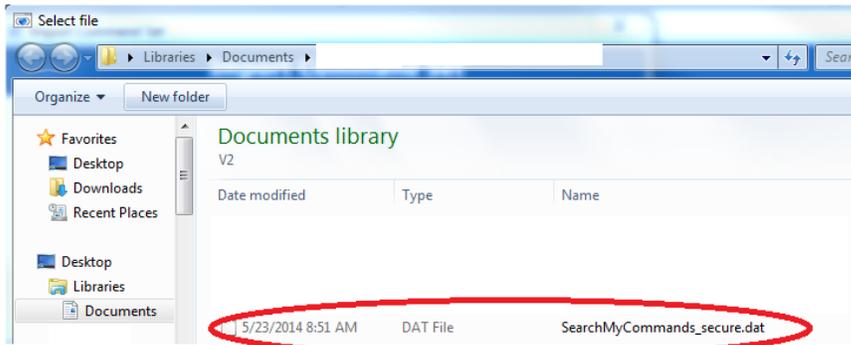

12. Click on Import

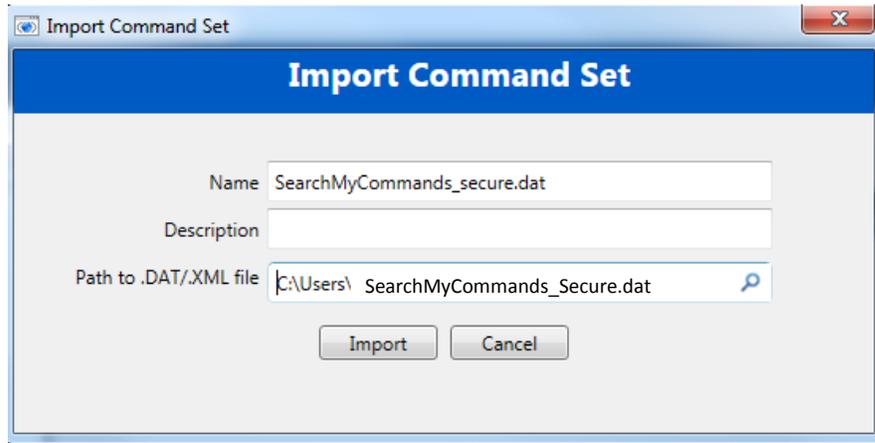

13. Click “OK” on the information window that appears letting you know the import has been scheduled.

- Wait for a few moments and click “Refresh” in the bottom right hand corner of the “Manage Command Sets” Tab

Command Set Name	Language	Owner Type	Owner Name	Description	Modified
SearchMyCommands_secure.dat	English (United States) [en-US]	Group	!Nuance Internal	14.04.09	6/24/2014 6:58 AM

Items: 5



- Repeat the above step until the command set appears in the window
- Repeat Steps 9 to 15 for each file to import.

Command Sets

- !Nuance Internal Site
- !Nuance Internal

[Click here to search...](#)

Command Set Name	Language	Owner Type	Owner Name	Description	Modified
.Calculators_1.11_secure.dat	English (United States) [en-US]	Group	!Nuance Internal		6/20/2014 7:28 AM
Commands.dat	English (United States) [en-US]	Group	!Nuance Internal		5/21/2014 12:31 PM
Conversations 2013 command examples.dat	English (United States) [en-US]	Group	!Nuance Internal		5/19/2014 12:55 PM
Nuance_PP-Content-Final-4-25-14.dat	English (United States) [en-US]	Group	!Nuance Internal	650+	6/12/2014 1:46 PM
SearchMyCommands_secure.dat	English (United States) [en-US]	Group	!Nuance Internal	14.04.09	6/24/2014 6:58 AM

Items: 5

